

ROGER MEDIA HUB CHEAT SHEET

CONNECTING TO ROOM AUDIO USING A MEDIA HUB

1. Turn on Media Hub

Disconnect the Media Hub from the power source.

Press and hold on/off button for at least 3 seconds until you get green flashing light on device.




2. Connect Media Hub to an Audio Source

2.1 In-built Hearing Augmentation system

Connect the Media Hub to the Wall Plate using the 3.5mm audio cable. Ensure that the sound is playing from the device with volume up.

2.2 Connect Directly to Audio Device

You can connect directly to an audio source by using the headphone jack of the computer/tablet, TV or phone. Ensure that the sound is playing with volume up.

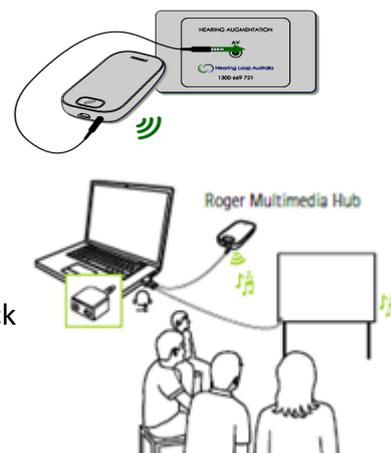


3. Connect Media Hub to a Hearing Aid

Press the link/connect button and hold within 5 cm to either the Roger Neckloop or Roger Receiver. A beep can be heard as the device connects.



4. Connect Media Hub to a Soundfield Speaker

Hold the Media Hub within 5 cm of the Roger Touchscreen. Press Connect button on the Touchscreen. Ensure that the Touchscreen is also connected to the Soundfield Speaker.

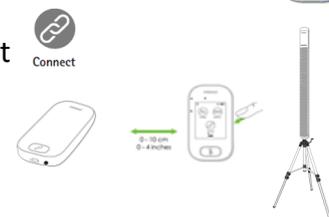

5. Reconnect the Media Hub to the power source

Troubleshooting

Tips

- Ensure that cables are connected to equipment with 3.5mm cables from the audio source ie. projector or TV.
- Check that the cable is connected at wall plate front and back.
- Ensure that the volume is up on all equipment that is being used.



	On/Off
	Charging
	Fully Charged
	Audio Broadcast
	Connecting
	Connection Successful
	Connection Failed
	Primary Device Not available (When paired with a touchscreen)



Hearing Loop Australia
Keeping your audience in the loop

www.hearingloop.com.au
info@hearingloop.com.au
1300 669 721