

PERSPEX BOX MOUNTING GUIDE

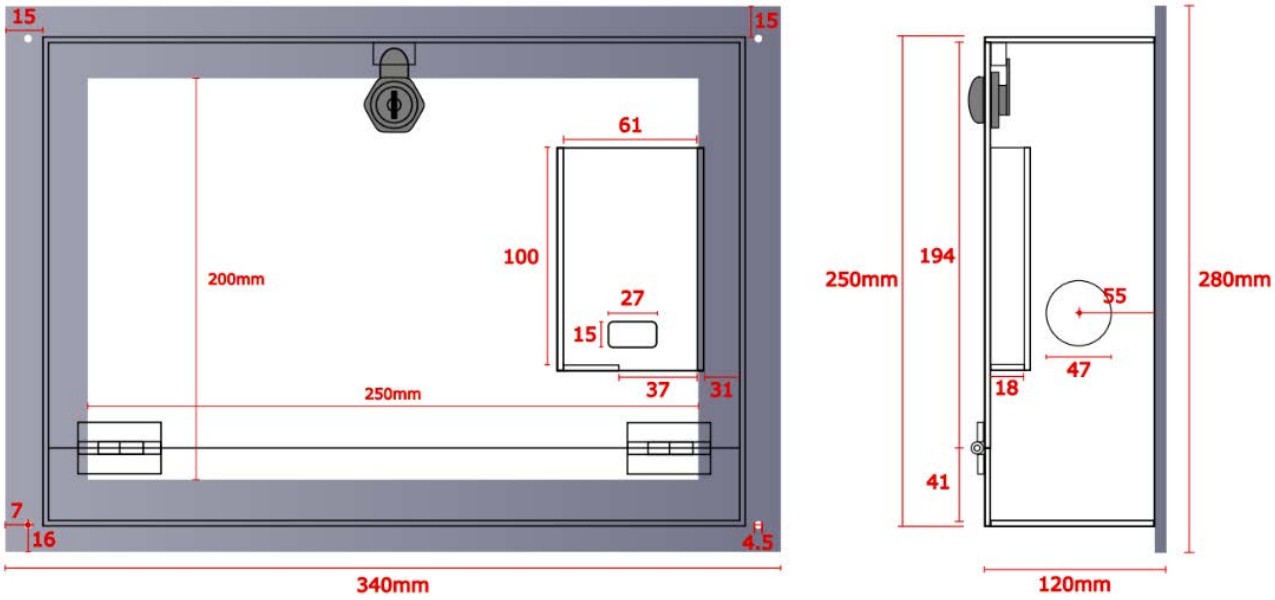

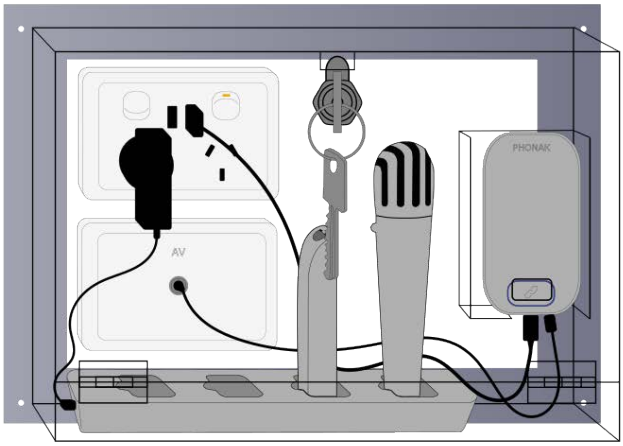

- Boxes are mounted using 4 screws (The type of screw will depend on the surface the box is being mounted to).
- Boxes NEED to be accessible as students may need to get their head within 10 cm to connect their hearing aids. We recommend a height of 1000mm AFFL.
- Boxes include an internal storage compartment for the Roger Media Hub, with front hole to access the connect button.
- Boxes have a 250mm x 200mm cutout in the back panel to enable clearance for various wall plate sizes that can be used for GPO or Audio.

ALL ELECTRICAL WORK SHOULD BE CARRIED OUT BY AN APPRPRIATELY QUALIFIED ELECTRICIAN

Hearing Loop Australia
Keeping your audience in the loop

www.hearingloop.com.au
info@hearingloop.com.au
1300 669 721

Connect & Communicate Pty Ltd t/as Hearing Loop Australia ABN 61 130 885 504